



Aligner Accessories

2017/18

Includes everything you and your patients need
to ensure your aligner treatments are a success

T: 0117 975 5533 **F:** 0117 975 5575
W: www.tocdental.com **E:** info@tocdental.com



Intensiv Ortho-Strip System



Safe and Precise

The Intensiv Ortho-Strip System combines a Kavo handpiece with Intensiv Ortho-Strips to provide the safest and most precise option available for interproximal reduction offering outstanding reliability. The Intensiv Ortho-Strip kit is an irrigated system reducing friction and heat thereby improving patient comfort. The comprehensive kit includes a Kavo contra angle handpiece with head, double sided Ortho-Strips, measuring gauge and cleaning block.

- Rapid and controlled reduction of enamel
- Safe and precise
- Outstanding reliability
- Fully autoclavable strips

Premium Kit

The reduction area may be approached from multiple angles as the handpiece head rotates 360 degrees and contouring is achieved using the flexible diamond strips.

Intensiv Ortho-Strip Kit Includes;

Kavo Contra Angle Handpiece Head, Set of 4 Double Sided Ortho-Strips, Set of Measuring Gauges and Diakleen Cleaning Block.

Code

KIT97-990-00



It is important to proceed systematically following a coarse to fine sequence in selecting grits. Unlike manual strips, the Ortho-Strip system provides rapid, controlled enamel reduction followed by polishing without unnecessary removal of healthy tooth substance.

Double Sided Ortho-Strips

- Flexible strips, diamond coated on both sides
- Proximal contouring, finishing and polishing of both adjacent teeth in a single procedure
- No injuries to patient (lips) or dentist (fingers)
- Available singly or in a set
- Fully autoclavable strips



Strips (Sold Individually)	Code
15µm Yellow - Pre-Polishing	97-990-19
25µm White - Finishing	97-990-18
40µm Red - Contouring	97-990-17
60µm Brown - Reduction	97-990-06
90µm Clear - Separating	97-990-16
Double Sided Set A (OS90/OS40/OS25/OS15)	97-990-12
Double Sided Set B (OS60/OS40/OS25/OS15)	97-990-11

Single Sided Ortho-Strips

- Flexible strips, diamond coated on one side
- For extra safety, these help preserve healthy tooth substance in the adjacent teeth
- No injuries to patient (lips) or dentist (fingers)
- Left and right versions for mesial and distal in the upper and lower jaw
- 4 different grades



Strips	Left	Right
15µm Yellow	97-990-25	97-990-29
25µm White	97-990-24	97-990-28
40µm Red	97-990-23	97-990-27
60µm Brown	97-990-22	97-990-26



Product	Code
Single Sided Set (L & R of each - OS15/OS25/OS40)	97-990-21

Incremental Thickness Gauges

Pack of six gauges for measuring interproximal spaces. Dimensions are 0.1mm, 0.2mm, 0.25mm, 0.3mm, 0.4mm, 0.5mm.

Code

IG-105-06



Diakleen

Rubber Block for rapid cleaning of clogged Intensiv Ortho-Strips and other diamond instruments.



Code

TOCACC-6000

Spacefile® Interproximal Reduction Kit

Versatility and Value

The SpaceFile Interproximal Reduction System is a complete system that includes the handpiece, the SpaceFile IPR files and accessories. The precision Swiss-made diamond-impregnated files can accomplish a wide range of tasks including interproximal reduction, contouring, polishing and finishing.

- Swiss-made, diamond-impregnated files
- Uses include breaking contact, reduction, contouring, fine trimming and finishing
- Designed to allow low temperature operation and prevent clogging
- Double-sided and single-sided
- Files, grip and inner tray are autoclavable
- Colour coded for easy identification
- Competitively priced



Spacefiles SWISS MADE



POLISHING

A final step in tooth surface restoration.



FINISHING

For finishing prior to final polish.



XX FINE

Pre-polish and finish reduced areas. Also use as initial strip for excessive crowding.



X FINE

Pre-polish stripped areas. Also use as initial strip for crowded cases.



FINE

Use to initiate reduction in crowded cases. Also use to contour after reduction.



MEDIUM

Use for reduction with some crowding or when some space exists.



COARSE

Primary working strip for IPR. Also use to trim and contour.



X COARSE

Use for maximum reduction at start of IPR.



SERRATED CONTACT POINT SAW

Use to break contact points on enamel, composite, amalgam, porcelain or gold.



REMOVABLE GRIP

Strips may also be used hand held with this removable grip.












Spacefile® Interproximal Reduction Kit Includes:

W&H Handpiece, Assorted pack of Spacefiles and Measuring Gauge.

Code

KIT04-500-10

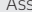Single Type Packs (Pk 8)		Double-Sided	Single-Sided (L)	Single-Sided (R)
	Polishing (DS .085mm, SS .070mm)	04-500-60	04-500-60L	04-500-60R
	Finishing (DS .110mm, SS .110mm)	04-500-30	04-500-30L	04-500-30R
	XX Fine (DS .140mm, SS .140mm)	04-500-15	04-500-15L	04-500-15R
	XFine (DS .160mm, SS .138mm)	04-500-25	04-500-25L	04-500-25R
	Fine (DS .270mm, SS .180mm)	04-500-40	04-500-40L	04-500-40R
	Medium (DS .330mm, SS .210mm)	04-500-150	04-500-150L	04-500-150R
	Coarse (DS .370mm, SS .235mm)	04-500-90	04-500-90L	04-500-90R
	X Coarse (DS .490mm, SS .290mm)	04-500-120	04-500-120L	04-500-120R
	Serrated Contact Point Saw	04-500-22	-	-
Assorted Packs		04-500-10	04-500-10L	04-500-10R

Single Type Packs include: Removable Grip and 8 files **Assorted Pack includes:** Removable Grip, 8 files (One of each thickness) and 1 Serrated Contact Point Saw

HandyStripper Kit

Compact Interproximal Stripping System


Specially designed for breaking contact points in tight interproximal spaces, the small HandyStripper delivers control and stability at your finger tips. Made of high grade autoclavable stainless steel with easy-to-replace diamond strips. Compact size so safer to use and variety of sizes makes it ideal for every interproximal stripping need.

Kit Includes: HandyStripper handle, 5 x Micro Saw 4mm, 10 x Single Sided 4mm Strips, 10 x Double Sided 4mm Strips

Product	Code
Handy Stripper Kit	IS-300-01
4mm Medium Single Sided Strip (Pk 10)	IS-300-02
4mm Medium Double Sided Strip (Pk 10)	IS-300-03
4mm Micro Saw Strips (Pk 10)	IS-300-04

Compact and lightweight

- Easy to control even in tight spaces
- Safe and controllable
- Durable made of stainless steel
- Autoclavable

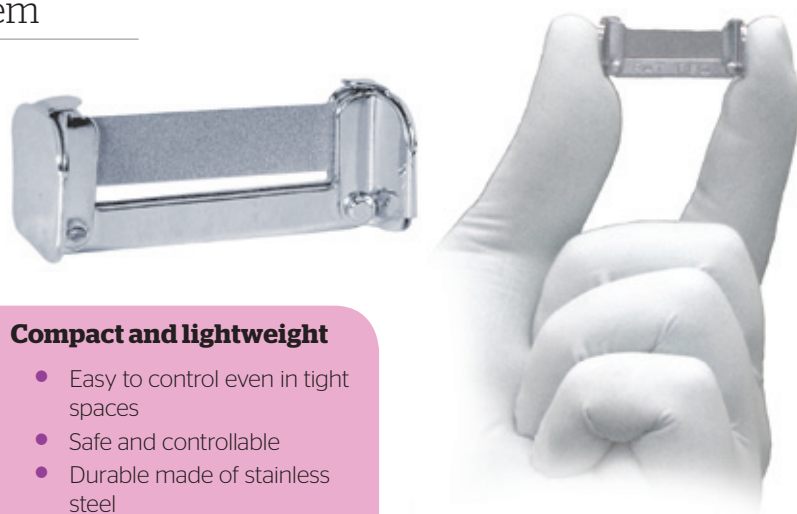

Intensiv ProxoStrip

Practical, fast, ergonomic



A unique, diamond-coated metal strip for final polishing of proximal surfaces. A metal strip with grip secure strip ends, diamond-coated on one side, with two working areas of different grit sizes.

- Sterilisable, reusable
- Finishing and polishing in one process
- Ergonomic grip for simple, safe use
- Short strip for improved finger support while stripping

Product	Code
Pack 6	TOC-PX4015-6
Pack 12	TOC-PX4015-12



Tooth Mousse

Use after interproximal stripping and after tooth whitening

- Reduces sensitivity
- Restores lost minerals
- Strengthens tooth enamel
- Neutralises acid pH



Assorted packs contain 2 x 90mg tubes of each flavour: mint, strawberry, tutti frutti, vanilla & melon.

Product (Pk10)	Code
Tooth Mousse, Assorted	TOC002032
Tooth Mousse, Mint	TOC002522

Anterior Air Rotor Stripping Diamond Burs

- No need to place a separator prior to using if teeth are parallel
- May be autoclaved, chemiclaved or dry heat sterilised - Don't immerse in cold sterilisation solution
- Friction grip



Product	Code
Head size 1.6mm, Head length 5mm	25-550-00

Air Rotor Stripping Diamond Burs

- May be autoclaved, chemiclaved or dry heat sterilised
- Do not immerse in cold sterilisation solution

Product	Code
Fine Grit	25-848-20
Medium Grit	25-848-10

Fine Grit



Medium Grit



Incremental Thickness Gauges

Pack of six gauges for measuring interproximal spaces. Dimensions are 0.1mm, 0.2mm, 0.25mm, 0.3mm, 0.4mm, 0.5mm.

Code
IG-105-06






Interproximal Diamond Strips

Perforated and flexible

- Perfect for interproximal reduction recontouring and finishing contact points
- Single sided
- Autoclavable
- 2 widths: Narrow (2.5mm) and Wide (3.75mm)
- 4 grades: coarse (90 micron), medium (45 micron), fine (30 micron), extra-fine (15 micron)

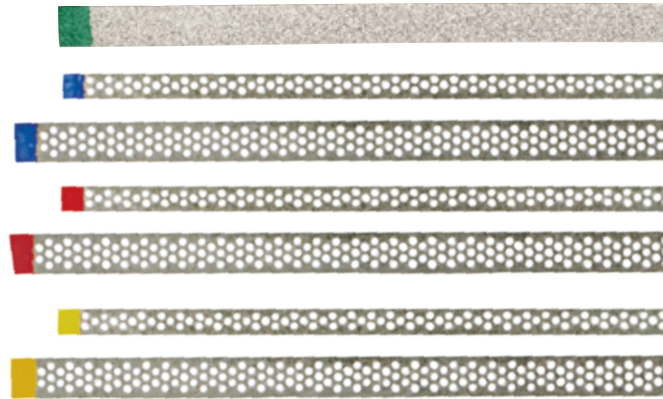

Interproximal Strip	Narrow	Wide
Extra Fine (Yellow Pk 10)	IS-600-03	IS-500-03
Fine (Red Pk 10)	IS-600-02	IS-500-02
Medium (Blue Pk 10)	IS-600-01	IS-500-01
Coarse (Green Pk 10)	-	IS-700-00
Narrow kit (Pk 15 - includes 5 each Extra Fine, Fine, Medium)	IS-600-00	-
Wide kit (Pk 15 - includes 5 each Extra Fine, Fine, Medium)	-	IS-500-00

Safe-Tipped Air Rotor Stripping Bur Kit

Sterilisable bur block containing 4 safe-tipped carbide burs, 2 medium particle size safe-tipped finishing diamond burs, 2 fine particle size safe-tipped finishing diamond burs, 2 anterior fine needle diamond burs, non-safe-tipped for anterior interproximal stripping.



Code

25-512-55

Interproximal Strip Tube

- Capable of sterilising 15 interproximal strips at a time
- Aluminium tube fits easily into your autoclave



Code

25-535-00

Interproximal Diamond Burs

It is recommended to begin with the medium diamond to quickly remove the required amount of enamel. Follow up with the fine diamond to smooth the tooth surface creating a precise, parallel end result.

- Consistent grit size
- Friction grip
- Head length 10mm
- Excellent diamond particle exposure

Heat sterilisable up to 340°F/170°C
Autoclavable up to 250°F/121°C



Product (Pk 5)	Code
Fine Grit (Red)	8859-31-010
Medium Grit (Blue)	859-31-010

Safe-Tipped Finishing Burs

- Diamond coated to help produce a natural finish
- Various particle sizes, all safe tipped
- May be autoclaved, chemically or dry heat sterilised - Don't immerse in cold sterilisation solution

Medium

Fine

Extra Fine



Product (Pk 3)	Code
Medium - Head 1.4mm, Length 10mm	25-540-55
Fine - Head 1.2mm, Length 8mm	25-530-55
Extra Fine - Head 1.2mm, Length 10mm	25-520-55

Diamond Discs

Smoothly Remove Tooth Mass

Flexible mesh design for both reduction and contouring. Smoothly remove tooth mass interproximally to create intra-arch space or open contact points between teeth. Great for use with aligner systems. Discs are available in 19mm, for when space is limited, and 22mm.



- Flexible for ease of use
- Disc diameter 19mm and 22mm
- For use with low speed handpieces
- Includes 2 diamond discs and 1 mandrel

Heat sterilisable up to 340°F/170°C / Autoclavable up to 250°F/121°C

Product (Pk 2)	Code
A. Double Sided Diamond Disc 19mm (0.17mm thick)	88617
A. Double Sided Diamond Disc 22mm (0.17mm thick)	88613
B. Double Sided See-Through Diamond Disc 22mm (0.20mm thick)	88655
C. Double Sided Mesh Diamond Disc 19mm (0.20mm thick)	88657
C. Double Sided Mesh Diamond Disc 22mm (0.20mm thick)	88656
D. Single Sided Diamond Disc 19mm (0.10mm thick)	88616
D. Single Sided Diamond Disc 22mm (0.10mm thick)	88614

Diamond Disc Safety Guards



- Use with Doriot style handpiece
- Disposable plastic option for complete visibility

Product	Code
19mm - Metal	88611
22mm - Metal	88612
22mm - Plastic (Pk 15)	88615

Diamond Disc Intro Kit



3 Disc Kit Includes:

1 x double-sided, 1 x single-sided,
1 x double-sided see-through
with mandrels and bur block

4 Disc Kit Includes:

1 x double-sided, 1 x single-sided,
1 x double-sided see-through,
1 x double-sided mesh with
mandrels and bur block

Product	Code
3 Disc 19mm Intro Kit	DD-KIT2
3 Disc 22mm Intro Kit	DD-KIT
4 Disc 22mm Intro Kit	DDKIT1

Heat sterilisable up to
340°F/170°C
Autoclavable up to
250°F/121°C

Optragate® Extra Soft Lip and Cheek Retractor

- Single use
- Flexible and elastic
- Latex Free
- Economical

Product (Pk 80)	Code
Assorted	TOC577275
Junior	TOC591451
Regular	TOC590850
Small	TOC590851



Lip Ring Retractor

- Protects lips during interproximal stripping
- Great for bonding and intraoral photography
- Made of high quality medical grade plastic

Code (Pk 2)

LR-400-02



Cold Sterilisable

Extraoral Cheek Retractors

- Provides a full view and easy access
- Strong one piece design

Sterilisable up to 370°F/188°C

Product	Code
Pedo	LC-100-01
Adult	LC-100-02



Extraoral Cheek Retractors

- Provides a full view and easy access
- Cold Sterilisable



Product (Pk 2)	Code
Pedo	LC-300-20
Teen	LC-300-21
Adult	LC-300-22

EXAFAST™

Fast Setting PVS Impression Material - Perfect for Clear Aligners

Exafast™ NDS Monophase



- Single phase technique with one viscosity
- Easy to syringe with enough compression to attain accurate impressions

Monophase Material: 2 x 48ml cartridges (96ml total), 6 mixing tips and instructions

Product	Code
Monophase Material	TOC137407
Mixing Tips II Large - Purple (Pk 48)	TOC800101

Exafast™ NDS Injection Type Wash



- Syringeable wash material for double phase impressions
- Thinner wash that is less thixotropic and more flowable
- Stays in place as pressure is applied to capture critical detail

Injection Type Material: 2 x 48ml cartridges (96ml total), 6 mixing tips and instructions

Product	Code
Injection Type Material	TOC137107
Mixing Tips II Small - Pink (Pk 48)	TOC800100
Intraoral Tips (Pk 100)	TOC800119

Exafast™ NDS Regular Type Wash



Most Popular

- Slightly thicker wash material used with double phase impressions
- Syringeable wash that is a little more thixotropic and less runny than injection material
- Stays in place until pressure is applied to capture critical detail

Regular Type Wash Material: 2 x 48ml cartridges (96ml total), 6 mixing tips and instructions

Product	Code
Regular Wash Material	TOC137007
Mixing Tips II Large - Aqua (Pk 48)	TOC900114
Intraoral Tips (Pk 100)	TOC800119

Exafast™ II NDS Registration

- Bite registration material that is easy to apply and will not slump or run
- Sets quickly in only 45 seconds and is easy to trim if necessary
- Ideal cream-like consistency allows patients to occlude without deviation or distortion



Bite Registration Material:
2 x 48ml cartridges (96ml total)
12 mixing tips and instructions

Product	Code
Bite Registration Material	TOC890014
Mixing Tips II L - Aqua pk 48	TOC900114
Cartridge Dispenser Type II NDS	TOC152101

Fast formulation with outstanding accuracy



Exafast™ Putty

EXAFAST Putty is a vinyl polysiloxane material, consisting of a base and a catalyst that are easy to manipulate for accurate, void-free impressions. EXAFAST Putty produces crisp, extraordinarily smooth and precise impressions and can be used with either regular or injection wash material.

Easy to Mix & Use



- The highest compression tray material used with double phase impressions
- Extremely smooth and easy to hand mix, avoid contact with latex, oils and moisture when mixing
- 45 second working time, sets after only 2¼ minutes in the mouth



Includes: 1x jar putty base (500g), 1x jar putty catalyst (500g), scoops and instructions

Product	Code
Putty Standard Pack	TOC137307

TOC Impression Trays

Perfect impressions first time!

- Full depth tray ideal for orthodontics
- Unique finger pressure platform assists balanced placement
- Slots with inner support tags retain alginate securely
- Excellent rigidity for improved accuracy of impression
- Long, non-slip handle for a firm grip
- Smooth rounded edges for added patient comfort
- No tray adhesive required
- Full range of sizes to ensure the right fit

Provides great accuracy needed for Aligner systems



Size (Pk50)	Upper	Lower
1 - Red (Pedo)	RT-100-01	RT-000-01
2 - Orange (Extra Small)	RT-100-02	RT-000-02
3 - Blue (Small)	RT-100-03	RT-000-03
4 - Green (Medium)	RT-100-04	RT-000-04
5 - Yellow (Large)	RT-100-05	RT-000-05
6 - Purple (Extra Large)	RT-100-06	RT-000-06

Assorted Bag (Pk 40 U/L Sizes 2-5)

RT-600-00

Upper

Lower



TOC Intra-Oral Elastics

TOC's Flower Elastics are manufactured from the finest surgical grade natural latex and latex free rubber to ensure reliable force delivery.




Colour coded bags for easy ID

- Fresh and uncontaminated due to air-tight seal
- Packs of 20 patient packs, each containing 100 elastics
- Patient bag includes 100 lightly powdered elastics for easy use with no residue

Product (Pk 20)	Colour	Force	Size	Latex	Latex-free
Daffodil	Yellow	2oz	5/16"	E20-516	-
Rose	Red	3.5oz	3/16"	E35-316	E35-316LF
Iris	Blue	3.5oz	1/4"	E35-14	E35-14LF
Mint	Green	3.5oz	5/16"	E35-516	-
Fuschia	Pink	3.5oz	3/8"	E35-38	-
Marigold	Orange	4.5oz	1/4"	E45-14	E45-14LF
Thistle	Grey	4.5oz	3/16"	E45-316	-
Honeysuckle	Tan	4.5oz	5/16"	E45-516	-
Heather	Mauve	6oz	5/16"	E60-516	-

Discover The Difference

Established for over 25 years, TOC is a leading supplier of orthodontic products to hospitals, specialist orthodontic and general dental practices throughout the UK and Ireland.

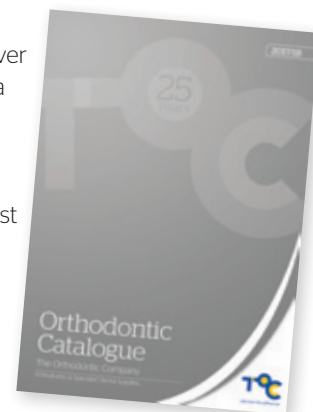

We offer a comprehensive range of products to suit all budgets and with one of the largest independent sales teams in the UK dedicated to orthodontics we are able to work closely with you to help understand your needs and assist in product choice.

If you would like a representative to visit please contact us.

TOC offers a complete range of Orthodontic Products - call us on **0117 975 5533** to receive your free catalogue or view it on our website **www.tocdental.com**

Freddy™ Ceramic Buttons & Hooks

The Freddy ceramic button is the perfect clear utility button for all aesthetic appliances.

- For labial and lingual use
- Mesh base creates a strong bond
- Bonds well with most common adhesives
- Anatomy and size of base fits all teeth
- Enough space around tie wing for all elastics
- Maintains the transparent look for aligners as anchorage for elastics

The Freddy hook has a smooth low profile and ellipse shape. It has the same characteristics as the Freddy button and matches perfectly with Invisalign™ and other aligner systems.



Product (Pk10)	Code	Product (Pk10)	Code
Bondable Ceramic Button	JB-334-00	Bondable Ceramic Hook	JB-335-00

Bondable Composite Buttons

- Aesthetic and stain resistant
- Contoured mechanical base for optimal bonding
- Ideal for use with aligners



Code (Pk 10)
JB-234-00

Bondable Composite Eyelets

- Aesthetic and stain resistant
- Contoured mechanical base for optimal bonding
- Funnelled entry for easy wire insertion



Code (Pk 10)
JB-235-00



Light cure your own aesthetic attachments - in seconds!




FILL the selected Mini-Mold Tip with light cure resin



PLACE the Mini-Mold Tip on desired tooth



LIGHT CURE the Mini-Mold Tip for 10 seconds

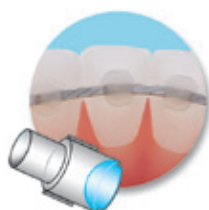

Wire Bonder Tips

Ideal for making direct bond lingual retainers. Comes in large and small sizes.



Posterior Molar Pad Tips

Designed to open the bite from the posterior teeth.



Bite-Ramp Tips

Creates a bite-opening ramp on the lingual of the upper anterior teeth.



Button Tips

Creates a functional low profile button quickly and easily.



Tube Tips

Designed to create a temporary tube.



Habit-Control Tips

The small, conical spike ends thumb sucking and tongue thrusting.



Bracket Tips

Creates edgewise brackets, ideal for minor tooth adjustments.



Starter Kit Includes

27 Tips, 3 of each style, 2 tip handles and detailed step-by-step instructions.

Type of Tips	Code
Bite Ramp 5mm (Pk 8)	51389
Posterior Molar pad (Pk 10)	51390
Tube (Pk 10)	51391
Bite Ramp 3mm (Pk 8)	51392
Wire Bonder Large 4.5mm (Pk 10)	51393
Button (Pk 10)	51394
Habit Control (Pk 10)	51395
Bracket (Pk 10)	51396
Wire Bonder Small 3mm (Pk 10)	51397
Spare Handle	FN-PLAS-133
Mini Mold Starter Kit	51399

NeoBond™ Light Cure Adhesive Kit

If your treatment plan requires the use of bondable attachments there is no better adhesive to use than NeoBond. The top selling adhesive from Dentsply is designed to work with LED lights as well as traditional curing lights. NeoBond works both chemically and mechanically for a superior bond.

Intro Kit



Product	Code
NeoBond Intro Kit	52-800-01
NeoBond Full Syringe Kit	52-800-00
4gm paste syringe	52-800-10

Ortho FlexTech™ Gold Chain Retainer Wire

Save on chairtime and lab fees. Clinically proven, Ortho FlexTech naturally forms to the lingual arch curvature making it easy to fit in the chair. The chain lays flat against the teeth offering improved patient comfort and lasting aesthetics due to its 14 karat gold composition.

An etched version is available for improved adhesion.

- Stabilises orthodontic corrections
- Prevents rotational relapses
- Low failure rate due to "flex linkage"



Product	Code
Gold Chain, 30" Spool	TOC-OFT
Gold Chain Etched, 30" Spool	TOC-OFT-E

Stainless Steel Chain Retainer Wire

An economical version of the gold chain with the same great patented design and ease of use. It offers a 30% increase in tensile strength and 300% increase in torsional strength.

Product	Code
Stainless Steel Chain, 30" Spool	TOC-OFTSS



FlowTain™ Flowable Light Cure Composite



FlowTain's low viscosity makes it ideal for bonding lingual retainers, splinting materials and indirect bonding of brackets. It is very smooth when cured which will reduce wear and enhance patient comfort. Save money and order this as a kit with OrthoFlexTech Gold or Stainless Steel Chain.

Product	Code
FlowTain 1.5gm Syringe + 10 tips	TOC-FIT
FlowTain + Gold Chain	TOC-OFTKIT
FlowTain + Gold Chain Etched	TOC-OFT-EKIT
FlowTain + Stainless Steel Chain	TOC-OFTSSKIT

Multi-Stranded Retainer Wire

Flattened Multi-Strand Wire in convenient 6" lengths. This is bonded to the lingual surface of the teeth as a permanent fixed retainer once orthodontic treatment has finished.

- No compliance issues
- Long term stability
- Improved patient comfort



Code (Pk 10)
FF-000-06

Pure Plastics™

Quick - Easy - Affordable



Opal™

Durable Retainer Material

- Excellent durability and strength
- First choice for retainers.

Product (Pk 100)	Code
1mm / .040" Circle - 120mm	O-40-120C
1mm / .040" Circle - 125mm	O-40-125C
1mm / .040" Square	O-40-127S

Diamond™

Whitening Tray Material

Offer your patients tooth whitening after treatment

- Easy to work with
- Exceptional performance



Product (Pk 25)	Code
1mm / .040" Circle	D-40-125R
1mm / .040" Square	D-40-127S
1.5mm / .060" Circle	D-60-120R
1.5mm / .060" Circle	D-60-125R
1.5mm / .060" Square	D-60-133S
2mm / .080" Square	D-80-133S

Ice™

Retainer Material

- Forms clear for great aesthetics
- Perfect for invisible retainers, minor tooth movement and temporary bridges

Product (Pk 20)	Code
0.75mm / .030" Square	I-30-127S
0.75mm / .030" Circle	I-30-125C
1mm / .040" Circle	I-40-120C
1mm / .040" Circle	I-40-125C
1mm / .040" Square	I-40-127S
1.5mm / .060" Square	I-60-127S

Easy Vac

Vacuum Forming Machine

Make your own whitening trays, retainers and gumshields with this quick and easy to use, stylish 'suck-down' machine.

- Compact size
- Easy to clean
- Safe and convenient controls
- Uses square material



Code
TOCEV2

Effer-ves®

Tablets for cleaning removable dental appliances

Provide your patients with a good cleaning regime from the start with the help of effer-ves. Used on a regular basis, effer-ves kills bacteria that can build up on aligners and other dental appliances as well as removing common stains caused by tea, coffee and some foods.

**Suitable
for use with
all aligners**



- Suitable for Clear Aligners, Clear Retainers, Wire Retainers, Mouthguards, TMJ Appliances, Bruxing Appliances and Nightguards
- Easy to use
- Safe to use on acrylic and wire combinations within orthodontic appliances

**Independently
tested
by Invisalign™**

Product	Code
Patient Box (32 tablets/box)	TOC-ACT
Poster	TOC-ACTP
Leaflet Holder	TOC-ACTLH
Patient Leaflets	TOC-ACTPL

Outie Tool

Easily and safely removes aligners and retainers

In order for an aligner to work correctly they need to have a tight fit, this sometimes causes a problem when the patient needs to take them out for eating and cleaning. Patients have reported sore fingers, broken finger nails and the use of unorthodox tools that could lead to damaging the aligners or worse, a tooth.

The special design of the Outie Tool allows the patient to remove the appliance safely and hygienically. The lever end is used to push the aligner up and off while the hook end is to pull the upper aligner down and off.

"My Outie Tool is SO easy to use. In just a couple of seconds, I dislodge my aligners high enough to be able to grab them and pull them out. It couldn't possibly be easier than this!"



Outie Tool
Retainer and Aligner Remover

Code (Pk 7)

PO-000-07

Easy Lift

A handy tool to have around the surgery. Unlike other versions of an aligner remover this one is made from high quality steel. If your patient has a gagging reflex this can help you pull off the aligners without having to put your fingers in their mouth before an examination



Code

TOC-EL

Retainer Retrievers

Retainer Retrievers help your patients remove their appliances from the buccal side, both cleanly and efficiently. They are a handy size and can be stored in an appliance case or on a key chain with the hole provided. Hygienically packaged in individual bags.



Product (Pk10)	Code
Tropical Blue	97-000-00
Tropical Pink	97-000-01

Chewies

Aligner Chewies are made of a soft plastic and are the same size as a small cotton roll.

Their purpose is to help the patient seat their aligner properly. Once the patient has put their appliance into the mouth they place the roll between the back teeth and bite down. This will help to push the aligner into the correct position giving the tight fit which is necessary for the treatment to work.



Product (Pk 20)	Code
Green - Mint	17-111-15
White - Unscented	17-111-16
Purple - Grape	17-111-17
Pink - Bubblegum	17-111-18
Yellow - Pineapple	17-111-19

Supplied in 10 individual bags of 2

ALIGNER CARE KITS

Essential Care for Aligner Wearers

Ensure all your patients take away one of these handy care kits with all the tools needed to look after their aligners during treatment. The durable bag is large enough to include additional items if needed.

Why not personalise these kits with your practice details - available at no extra charge when you order a minimum of 24 kits.

Please note: Aligner case and Outie Tool colours may vary. Bag may vary slightly to example shown

Kit contains:

1 x box of Effer-ves®
& Patient Leaflet
Cleaning tablets

1 x Aligner case
For hygienic and safe storage of aligners

1 x Pack of Chewies

To help seat aligners correctly and comfortably in the mouth

1 x Outie tool

For easy removal of aligners and retainers from the mouth

NEW



Product	Code
Aligner Care Kit	PA-100-00

Appliance Cases

Slim design ideal for aligners, whitening trays and retainers. Try our pack of assorted colours so you have a colour to suit everyone!

- Strong & compact
- 1" Depth
- Vented
- Economical
- Positive lock hinged lid

GREAT VALUE



Product (Pk 12)	Code
Pink	TOC-RCPK
Blue	TOC-RCB
Coral	TOC-RCC
Green	TOC-RCG
Yellow	TOC-RCY
Orange	TOC-RCO
Purple	TOC-RCP
Red (Pk 10)	TOC-RCR
White	TOC-RCW
Assorted	TOC-RCA
Assorted Sparkle	TOC-RCS

Retainer Cases

Strong, high-impact plastic in your choice of classic, glow in the dark and sparkle finishes with write-on I.D. area

- Appliance key holder style
- Durable long life hinge
- 1 1/4" depth fits most appliances
- Vented holes prevent odours

BEST SELLER



Product (pk 25)	Code
Blue	350-000
Red	350-001
Pink	350-002
Orange	350-003
Yellow	350-004
Green	350-005
Black	350-006
Clear	350-007
Aqua	350-008
Purple	350-009
Assorted	350-010
Assorted Glow	350-020
Assorted Sparkle	350-030



Deluxe Appliance Case with Mirror

When your patients want something a little more stylish to pop in their handbag or suitcase. A smart sturdy appliance case with a handy mirror within the lid.

- Perfect for aligners and retainers
- Smooth rounded compact design
- Secure snap-tight closure
- Integrated metal reinforced hinge
- Durable plastic with high gloss polished finish
- 3" x 1" with 3/4" inside depth



Product (Pk 12)	Code
Black & White	45204
Assorted Colours	45205